


1

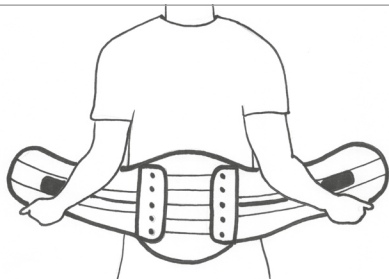
2

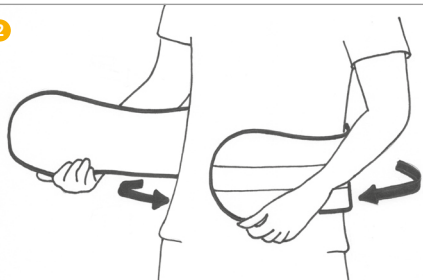3

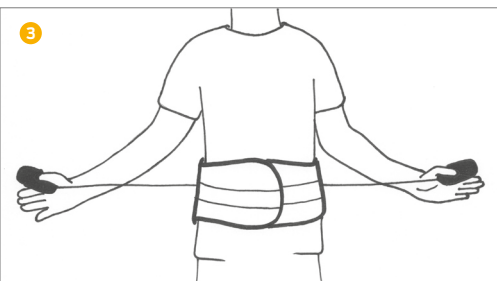4

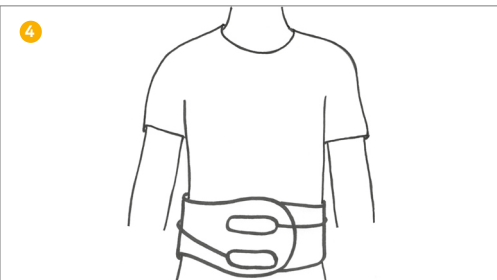DK

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Intended use:
Lumbar orthosis for the support, strain relief and stabilisation of the lumbar spine.

Indications:
Lumbago
Lumboischialgia
Herniated disc
Spondylolisthesis
Stable lumbar vertebra fractures
Postoperative stabilisation after intervertebral disc surgery

Contraindications:
Currently unknown.

Application instruction:
Open the front fastener of the orthosis and release the pull cords. Place the orthosis around the waist (1). Position the orthosis so that the back section lies centred in the lower back area. If necessary, the two rigid side pieces can be repositioned for a better fit using the Velcro fastener. When repositioning, make sure that the side pieces rest without pressure. Close the two sides of the front fastener evenly over one another using the Velcro fastener (2). Grasp the two lateral pull cords and pull them forward with an even traction (3). Once the desired stabilisation has been achieved, fasten the pull cords using the Velcro fasteners (4). To take off the product, first release the lateral pull cords and then open the front fastening.

- Application note / Caution:**
- Users are advised not to wear the product directly on the skin, but for example over a shirt.
 - Only change the settings of the product after consulting your attending physician.
 - Do not use the product on open wounds. In the following cases you should consult your doctor before use: Skin diseases, lymphatic drainage disorders, sensory disorders, circulatory disorders, unclear soft tissue swelling.
 - This product is intended for use on a single patient only.
 - If necessary, the help of a second person can be beneficial for the correct application of the product.
 - Only by choosing the correct size and variant, as well as the correct application, an optimal function can be achieved.
 - The product should no longer be used if it is damaged or heavily worn, or if existing joints and fasteners no longer function properly.
 - Make sure that existing Velcro fasteners do not cause damage to clothes or other fabric parts.
 - Do not put on the product too tightly, as this may affect the blood circulation.
 - If symptoms persist or worsen, please consult your medical doctor.
 - Please inform the manufacturer and the competent authority in your country about serious incidents with this product.
 - If you are unclear or have questions about the product and its handling, please consult your doctor or qualified personnel.
 - Do not make any changes to the product.
 - Simultaneous use with other products only after consultation with your doctor.

Care instructions:
The textile parts of the product can be cleaned by hand with water and a mild detergent. The plastic and metal parts of the product can be cleaned with a damp cloth. Close existing Velcro fasteners before washing to avoid damaging other laundry items. Allow to air dry. Avoid direct heat sources and do not machine dry or iron. Make sure the product is completely dry before using it again. Do not dry clean and do not use bleach, solvents, or fabric softeners. Fatty or acidic agents, ointments or lotions can damage the material.

Material composition:
Polyester (PES), PUR foam, ABS plastic, Aluminium, Polyoxymethylene (POM), Polyamide (PA), PVC

Storage instructions:
Store the product in a dry place and protect it from heat, direct sunlight and moisture.

Disposal note:
Disposal via household waste is possible.

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